

Unique and practical solutions in ostomy and woundcare



Ostomy care

Wound care

Accessories

CE

1 • Skin wafer

The skin friendly and very flexible skin wafer is made of 100% hydrocolloid. It provides a secure fit, protects the skin and absorbs moisture efficiently. There is excellent balance between secure adhesion and comfortable removal. The skin wafer is available in pre-cut sizes or can easily be cut to the desired size using the cut-to fit version.




2 • Tissue

The very soft, breathable and water resistant tissue is comfortable against the skin and prevents perspiration forming underneath the pouch. It dries very quickly after swimming or showering. The soft film layer has an excellent oxygen barrier, totally preventing odour transmission.



3 • Filter

Our exceptionally good activated carbon filter provides a very long odour neutralising pathway and can cleanse large amounts of gas. This means the gas will be odour free when it leaves the filter. This filter also prevents ballooning of the ostomy pouch because it has good air flow; however, leakage through the filter is prevented by the protective inner layers of the pouch and the long odour neutralising pathway within the filter.

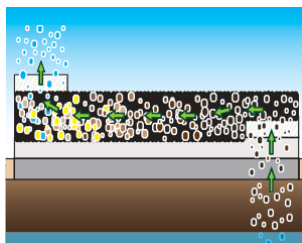

4 • Viewing window

The ostomy pouch includes a viewing window which makes it possible to inspect the stoma as well as the output from the stoma. It is also very helpful when applying the pouch to make sure it is located correctly over the stoma.



5 • Integrated vent

The integrated vent is a new and novel way to vent the pouch allowing air out or in, for instance in the case of a vacuum occurring. It can also be used for more effective cleansing by rinsing the pouch via the vent opening.



6 • Clamp or integrated velcro type closure

Ileostomy pouches can be securely closed by using a clamp or soft tie, both of which are provided. Our innovative range of ileostomy products also includes a pouch with an integrated velcro type closure to use instead of a clamp or ties, making the opening and closing of the pouch a much easier task. This strong integrated velcro type system can conveniently be closed by rolling up the outlet and securing the closure via the hook and loop velcro type pads attached to the film. Opening the pouch is equally as easy by disconnecting the two velcro type pads and unrolling the outlet. Additionally, at the base of the outlet is a plastic strip that acts as an opening aid, making it much easier to hold the pouch open when emptying.




7• 2-piece ring system

The reliable, secure, and solid ring system is very flexible which provides high comfort to the user. The pouches can be easily attached to the base plate. The floating ring ensures a very strong and flexible connection between the pouch and the base plate.



8• Non-return valve

Due to the non-return valve preventing the urine from flowing back towards the stoma, the risk of leakage is minimal even when lying down and comfort to the user is increased. The urine bag has an innerlayer: when the bag is upside down, the urine cannot flow back to the stoma.

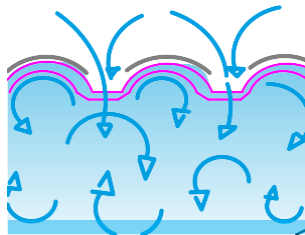

9• Urine tap

The secure and user friendly urine tap can easily be connected to a night drainage bag or leg bag without using an extension connector. This suits to the most standard connectors.




10• Ultra Light

The ZenSiv Ultra Light skin wafers are the answer to leakage problems. The skin wafer is extra thin and has very flexible and comfortable touch and feel. The Ultra Light wafer has ribs for extra secure adhesion on the sides, the ribs avoid curling of the edges.



11• Foam layer

Some of the ZenSiv hydrocolloid skin wafers are finished with a very soft foam backing. This soft and smooth foam layer makes the wafer and pouch very comfortable to wear.

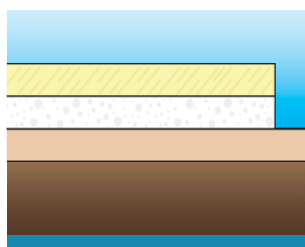

12• Irrigation of the ileostomy pouch

The integrated vent in our drainable pouches allows the pouch to be irrigated during wearing. A hygienic way to clean the stoma and pouch. More hygiene and comfort for the user due to this unique solution.



Our adhesive strength

Having a stoma can lead people to feel insecure about themselves. They are concerned about the pouch staying in place or leakage. Feeling secure about your bag is something every ostomist desires. ZenSiv has developed an adhesive that is adapted to that: the adhesive strength of our base plate has been deliberately strengthened.

Because of the high initial attack the wafer sticks to the skin in one go. And all without difficult instructions or having to apply supplementary aids. Security above all. So choose the adhesive strength of ZenSiv!

Security:

- Comfortable skin wafer with excellent adhesion: for secure and guaranteed adhesion
- Extra secure adhesion, even on the edges thanks to the ribs



Our filter and air vent

Ballooning or vacuum formation of the pouch and escaping of unpleasant odours, all feared problems for people with a stoma. Research is being carried out worldwide into these issues. ZenSiv offers a unique solution to ostomists. Our unique, integrated vent makes it possible to allow air into the bag or be released discretely and hygienically. Nobody will notice a thing!

The durable filter in the ZenSiv range prevents the dispersion of unwanted odours. You can wear your bag without concerns, without air problems.

That is what every ostomist desires.

Control:

- Unique integrated vent: you control the air
- Durable filter: control of odours



Our skin comfort

Skin-friendliness is a commonly used term when talking about ostomy products. But only skin-friendliness is not good enough for us. We go one step further and talk about skin comfort, since no ostomist wants to have red or irritated skin. Why compromise on adhesion?

The ZenSiv skin wafer offers natural skin comfort. Effortlessly adapted to you and comfortable as it should be. The material has been well tested and some ostomists even forgot they were wearing it. That is skin comfort according to ZenSiv!

Comfort:

- The comparative of skin-friendliness: skin comfort
- Extremely flexible and comfortable



Our products do not contain any animal ingredients out of considerations for humans and animals.



OxMed is campaigning against leakage. Our hydrocolloid is specifically designed to prevent leakage. The direct and strong adhesive strength prevents leakage and provides a more secure protection right from the beginning.

A complete solution for control, security and comfort? This is possible with the unique solutions from ZenSiv!

ZENSIV
OSTOMY
products

1-piece and 2-piece colostomy



1-piece and 2-piece ileostomy



1-piece and 2-piece urostomy



ZenZipp All-In-One Wound Management System

The traditional care for large, difficult to manage (high-output) wounds is time consuming and challenging for healthcare professionals and often very uncomfortable for patients. The ZenZipp All-In-One wound management system was developed in close cooperation with specialists, nurses and patients to improve this situation:
A unique, complete and practical wound management system with countless advantages.

ZenZipp All-In-One Wound Management System consists of a line of flexible and user-friendly wound managers in different sizes, suitable for every type of wound and every body type. The product is fitted with a flexible and skin-friendly hydrocolloid sticky layer for optimal adhesion and skin protection.

ZenSetiv hydrocolloid paste and hydrocolloid strips can be used to fill skin irregularities. This creates an optimal seal for the wound, preventing leakage and unpleasant smells. ZenZipp All-In-One wound managers can stay in place for several days.

Thanks to the transparent window in the sturdy collection bag, the wound can be inspected properly at any time without interfering with the healing process. An easily opened and closed double zipper over the entire length of the bag allows quick and comfortable wound care. Finally, the taps allow effective draining. These taps allow problem-free connection with catheter bags, a vacuum pump or other collection bag.

The advantages of ZenZipp All-In-One Wound Management System

• COMFORTABLE CARE

The flexible wound managers are leak-free and can stay in place for several days: Greater comfort for the patient and significant savings of both money and time.

• TRANSPARENT ACCESS PANEL

The large, transparent access panel with integrated double zipper allows full access to the wound for observation and treatment, without removal of the wound manager.

• OPTIMAL ADHESION

The soft, breathable hydrocolloid adhesive ensures optimal adhesion and flexibility, high comfort levels and prevents skin reactions.

• AIR CUSHION

Through the integrated air vent it is possible to insert a large amount of air to prevent the film from touching the wound area. The wound healing process will therefore be unharmed and leads to an increase of comfort for the patient.

• VELCRO STRIPS


When opening the pouch on the four larger sizes, the edges of the opening can be fixed on the adjacent Velcro strips so the bag stays open without the need to hold it, allowing easy access with both hands.

• DRAINAGE AND COLLECTION CAPACITY

The integrated collection bag provides good draining and collection capacity, specially designed for high-output wounds. For example a wound with a fistula, heavily infected wounds and/or a wound in combination with a stoma.

• DRAINAGE

Irrigation, drainage and catheterization are no problem thanks to the taps fitted to the bag. If a tap is not needed, it can be easily cut off or folded back: the integrated fold up closure provides a safe closure of the wound manager. It doesn't matter how the wound manager is placed: horizontally or vertically the content of the wound manager can easily be drained through one of the three taps. There is always at least one that is facing down for comfortable draining.

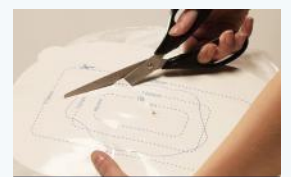

Optimal care
with a single ZenZipp wound manager;
practical, comfortable, time-saving
and cost-saving

ZEN
ZIPP.

- 1 unique solution
- 2 fit all wounds
- 3 taps for secure drainage

1, 2, 3 go!

A fit for all sizes

Available in different sizes, a suitable wound manager for every wound.
Individual sterile packaging, lightweight, easily cut to size.

Wound Managers

Maximum cut sizes
170x120mm; 220x170mm; 270x190mm and 390x240mm



Maximum cut size 110x75mm



Maximum cut size 110x110mm



Maximum cut size 50x30mm



ZenSetiv Hydrocolloid Paste

The ZenSetiv Hydrocolloid Paste is a filler and skin barrier that has been specially designed for use around stomas and fistulas.

The Hydrocolloid paste contains no alcohol, is sting-free and very skin friendly.

Advantages:

- No Alcohol
- Minimizes leakage between the skin and the skin wafer
- Extremely skin friendly



ZenSetiv Stomabelt

The ZenSetiv belt is made of soft and comfortable material, and is very discreet to wear. It ensures the position of the stoma pouch for extra security.

Advantages:

- Extra Security
- Soft material
- Washable



ZenSetiv Remover Spray

ZenSetiv Remover Spray is a specially formulated silicon based solution that ensures safe and pain free removal of all medical adhesive devices and skin residues.

The ZenSetiv adhesive remover is particularly suited for removing difficult residues which build up on the skin over time.

The Remover Spray is packed as a pumpspray and is very economical in use. The product does not contain alcohol and does not sting.

Advantages:

- Sting free
- Alcohol free
- Economical in use



ZenSetiv Protective Powder

ZenSetiv Protective Powder is used to absorb moisture from the skin around the stoma after cleaning the area when changing, in order to assist adhesion of a fresh stoma appliance.

Advantages:

- Can be applied easily
- Has a soothing effect on the painful skin
- Protects the skin





Affordable quality
products and care solutions

ZENsetiv

ZenSetiv Mouldable Hydrocolloid Rings

The ZenSetiv mouldable hydrocolloid rings are used to protect the skin around the stoma. The rings absorb moisture and provide a secure fit. By covering the skin around the stoma with a mouldable ring, the skin is very well protected against aggressive output. Skin irritations that are caused by aggressive output are avoided. The skin protection rings are mouldable and can be adapted to different types and sizes of stomas.

Advantages:

- Optimal skin protection
- Minimizes leakage between the skin and the skin wafer
- Mouldable: useable for different types and sizes of stomas

ZenSetiv Silicon Plasters

The ZenSetiv Silicon Plasters are the ideal solution for ostomists that irrigate their stoma or with a continent ileostomy. The soft silicon layer sticks to the skin immediately and protects the stoma. The foam layer absorbs mild secretion of liquid stool. The ZenSetiv Silicon Plasters are very comfortable and painless to remove. Also available as silicon plaster with split, which can be used to cover body openings around swabs, catheters, or gastro tubes for example.

Advantages:

- Ideal protection for ostomists that irrigate or with a continent ileostomy
- Fast initial tack
- Very skin friendly
- Absorbing foam layer
- Painless to remove





Quality assurance

OxMed International GmbH has quality registration with Lloyds Register Quality Assurance to ISO 13485:2003 and all products are CE marked. The products meet the highest quality standards from our choice of raw materials, through production to packing and despatch.

We only use the best materials and our products must regularly pass strict tests and quality control procedures. This is how we are able to ensure comfort and security in our product ranges.

Moreover, we can guarantee consistency which, in turn, assures confidence to the wearer.



Distributor:

